

SPECIFICATIONS

Model: LDBS3200

Measurement: Body weight

Top material: 6mm tempered safety glass

Bottom material: ABS plastic

Auto power ON/OFF: Yes

Low battery indicator: Yes

Overload indicator: Yes

Display: Red LED on black backlight

Display size: 80 x 32 mm

Capacity: 180 kg / 396 lb

Division: 0.1 kg / 0.2 lb

Power supply: 3 x 1.5V AAA batteries

Unit size: 300 x 265 x 24 mm

! ADVISE FOR USE & CARE

- 1.The platform will be slippery when it is wet. Keep it dry!
- 2.Keep standing still during the course of weighing.
- 3.DO NOT strike, shock or drop the scale.
- 4.Treat the scale with care as it is a precise instrument.
- 5.Clean the scale with a damp cloth and avoid water getting inside the scale.
DO NOT use chemical/abrasive cleaning agents.
- 6.Keep the scale in cool place with dry air.
- 7.Always keep your scale in horizontal position.
- 8.If failed to turn on the scale, check whether battery is installed, or battery power is low. Insert or replace new battery.
- 9.If there is error on display or can not turn off the scale for a long time, please take off the battery for about 3 seconds, then install it again to remove soft-fault.
If you can not resolve the problem, please contact your dealer for consultation.
- 10.Do not use it for commercial purpose.



L'Docile Enterprise GmbH Germany


Radlkoferstrasse 2 - 81373 - München, Germany

info@ledocile.com www.ledocile.com

L'DOCILE
HEALTHCARE PRODUCTS

Electronic Personal Scale


INSTRUCTIONS FOR USE

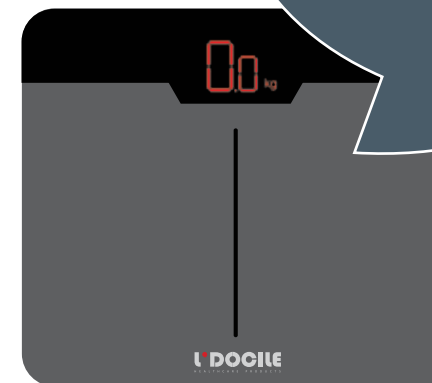

• The actual item may vary from the photographs/illustrations.


LDBS3200

BATTERY

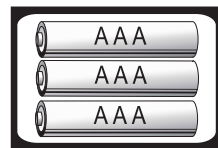


3 x 1.5V AAA batteries

OPERATION




Step 1: To make sure accuracy, place the scale on a hard, flat floor and avoid carpet or soft surface.

Remark: If the scale is one of our "slim" series (indicated on package), make sure to put it on hard surface such as marble, floor tile etc., and make sure no foreign bodies beneath the scale before weighing.

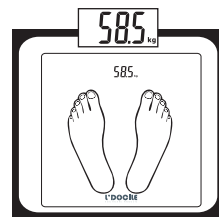


Step 2. Gently step on the scale, then the scale will be automatically turned on. Stand evenly on the scale without moving and wait until your weight shown on display is stable and locked.

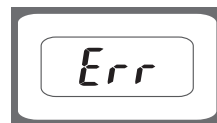

OVERLOAD INDICATION


Scale overload. (max. weight capacity is indicated on the scale).
Please step off to avoid damage.

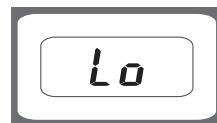

LOW BATTERY INDICATION


Battery power is low. Please install new battery.

TO REPLACE BATTERY

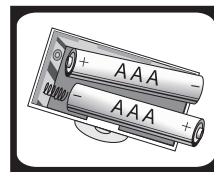

Open battery case cover on the back of scale
Take out the used batteries (with the aid of a cuspidal object, if needed)
Install new battery by putting one side of the battery and then pressing down the other side.
Please take note to the polarity.


AUTOMATIC SWITCH OFF

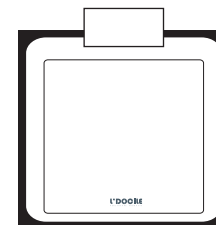

When you step off the scale, it will automatically switch off. Auto shut-off occurs if display shows the same weight reading for approximate 8 seconds.

TO CONVERT WEIGHT UNIT

If the scale is with both measuring system of Metric(kilograms) and Imperial (pounds, stones), you can convert weight unit as follow:

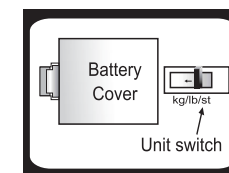

The unit switch is on the back of the scale
Find unit switch on the back of scale.
Push up and down the unit switch to convert scale unit between kg-lb-st.